


# Bicos Meniscal Repair Retractor

Designed by James Bicos, MD

*A popliteal retractor specifically designed for meniscal repair or access to the posterior knee*

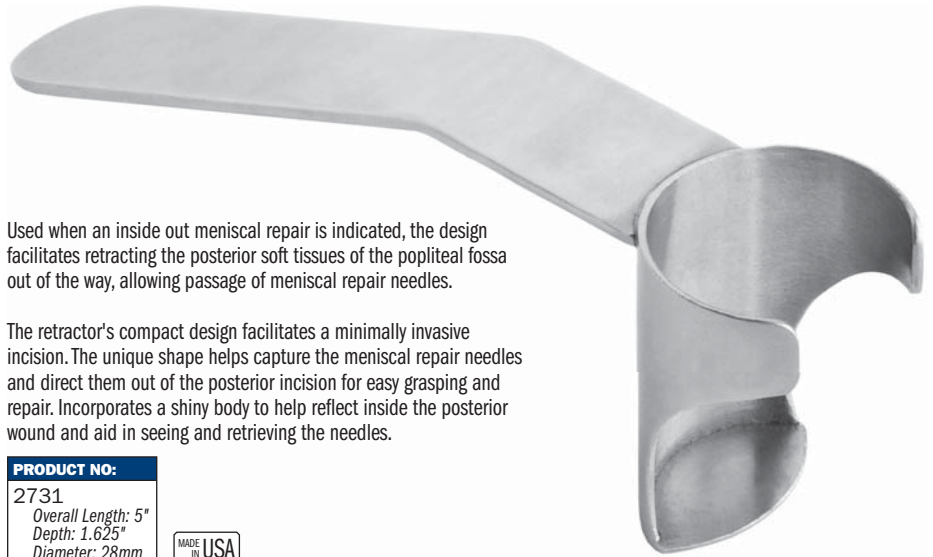

Used when an inside out meniscal repair is indicated, the design facilitates retracting the posterior soft tissues of the popliteal fossa out of the way, allowing passage of meniscal repair needles.

The retractor's compact design facilitates a minimally invasive incision. The unique shape helps capture the meniscal repair needles and direct them out of the posterior incision for easy grasping and repair. Incorporates a shiny body to help reflect inside the posterior wound and aid in seeing and retrieving the needles.

**PRODUCT NO:**

2731

Overall Length: 5"

Depth: 1.625"

Diameter: 28mm



FREE TRIAL ON MOST INSTRUMENTS

© 2012 Innomed, Inc.

**INNOMED** 



103 Estus Drive, Savannah, GA 31404  
www.innomed.net info@innomed.net

912.236.0000 Phone  
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 73  
Fax +41 41 740 67 74

1.800.548.2362