



# Curved Cement Osteotome

*Helps remove cement around the back of the tibia base, and is useful in the femoral notch during removal of a knee femoral component*

Designed to be inserted around the back of the tibia base to remove cement. The curve is congruent with most tibia bases. During revision knee surgery, can be used to help separate the prosthesis/bone or prosthesis/cement interface. The curve of the osteotome allows it to be used in the femoral notch of a femoral component. The osteotome is nitrate coated to help protect the implant surface.

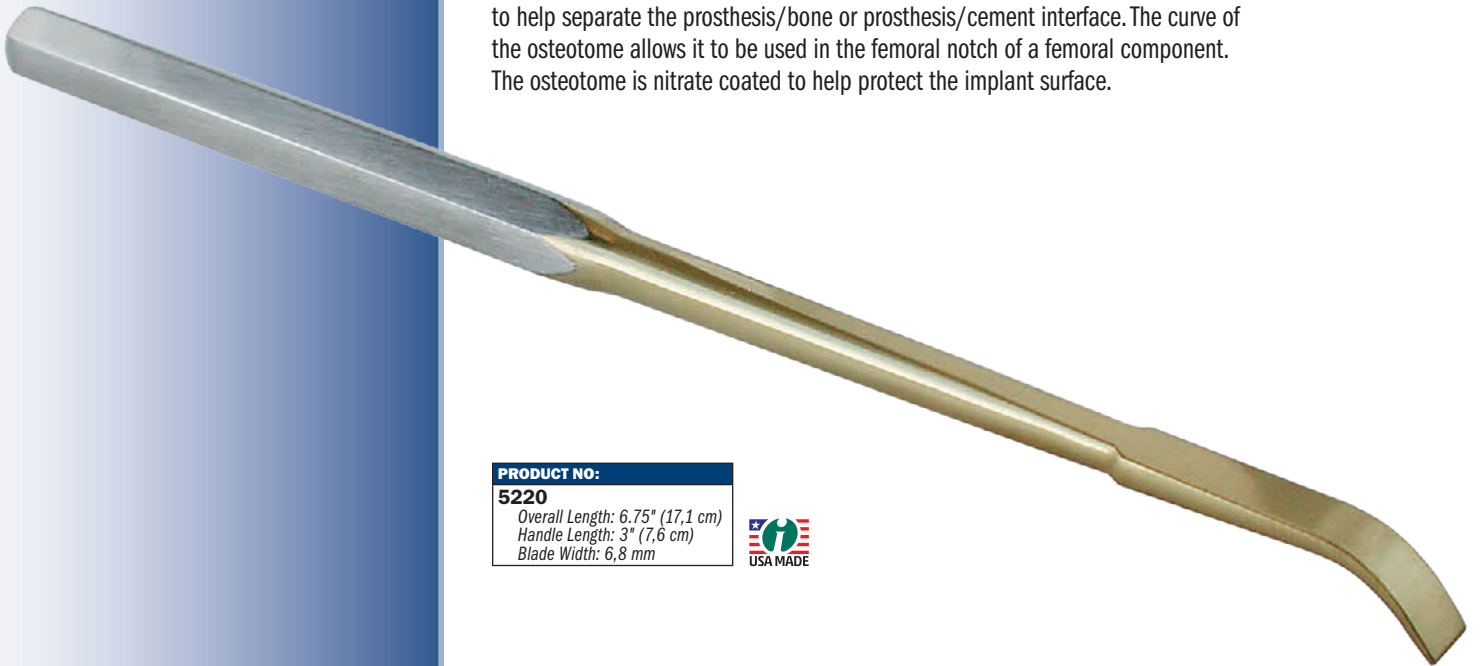

**PRODUCT NO:**  
**5220**  
Overall Length: 6.75" (17,1 cm)  
Handle Length: 3" (7,6 cm)  
Blade Width: 6,8 mm

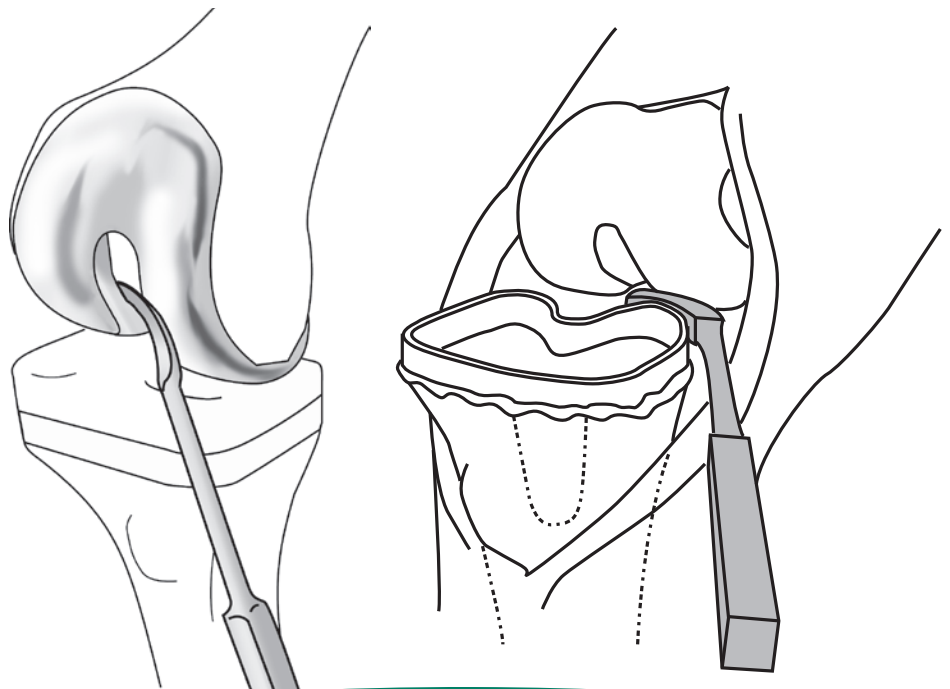

[www.innomed.net](http://www.innomed.net)

ISO 13485:2016

**FREE TRIAL ON MOST INSTRUMENTS**

# INNOMED



103 Estus Drive, Savannah, GA 31404  
info@innomed.net

912.236.0000 Phone  
912.236.7766 Fax

**Innomed-Europe** Tel. +41 41 740 67 74  
info@innomed-europe.com

**1.800.548.2362**