



Hagan Carpal Tunnel Release Sleeve



Designed by Hugh Hagan, MD

Designed to protect the surrounding anatomy while providing a sleeve within which to smoothly advance a beaver-style blade to divide and release the transverse carpal ligament

Designed for use in a mini-open, non-endoscopic approach, the sleeve isolates the blade, providing protection to the surrounding anatomy. The longer, bottom leading edge of the sleeve is inserted between the median nerve and the transverse carpal ligament, while the shorter, top leading edge provides lifting protection to the structures above the ligament. The blade is then advanced within the sleeve to complete the ligament release.

PRODUCT NO:

1150

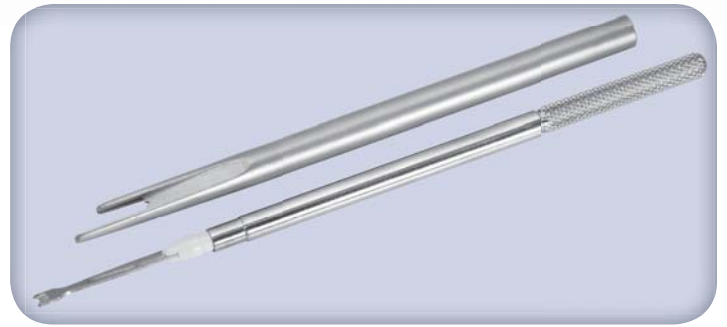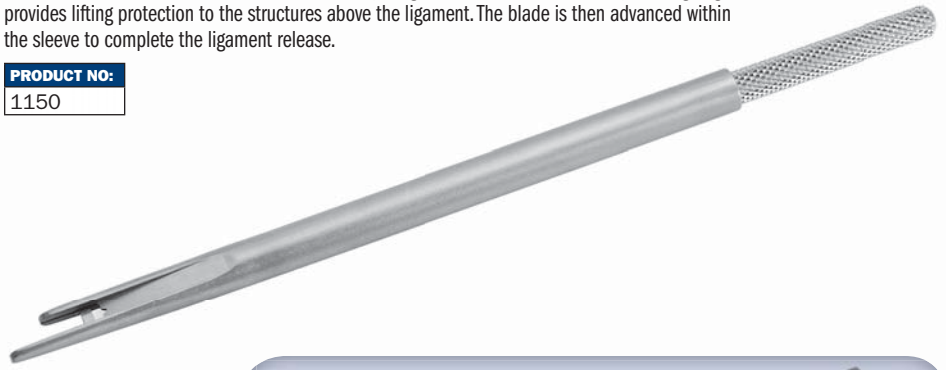

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax



INNOMED

1.800.548.2362