


# Kenerly Femoral Neck Cutting Guide

Designed by J. Lex Kenerly, III, MD

*Designed for use during the anterior approach for THA to help determine the femoral neck osteotomy location for the initial femoral neck osteotomy*

The guide is placed on the femoral neck and adjusted using the intraoperative C-arm image to visualize and compare to the pre-op templating, providing an excellent location for the initial femoral neck osteotomy.

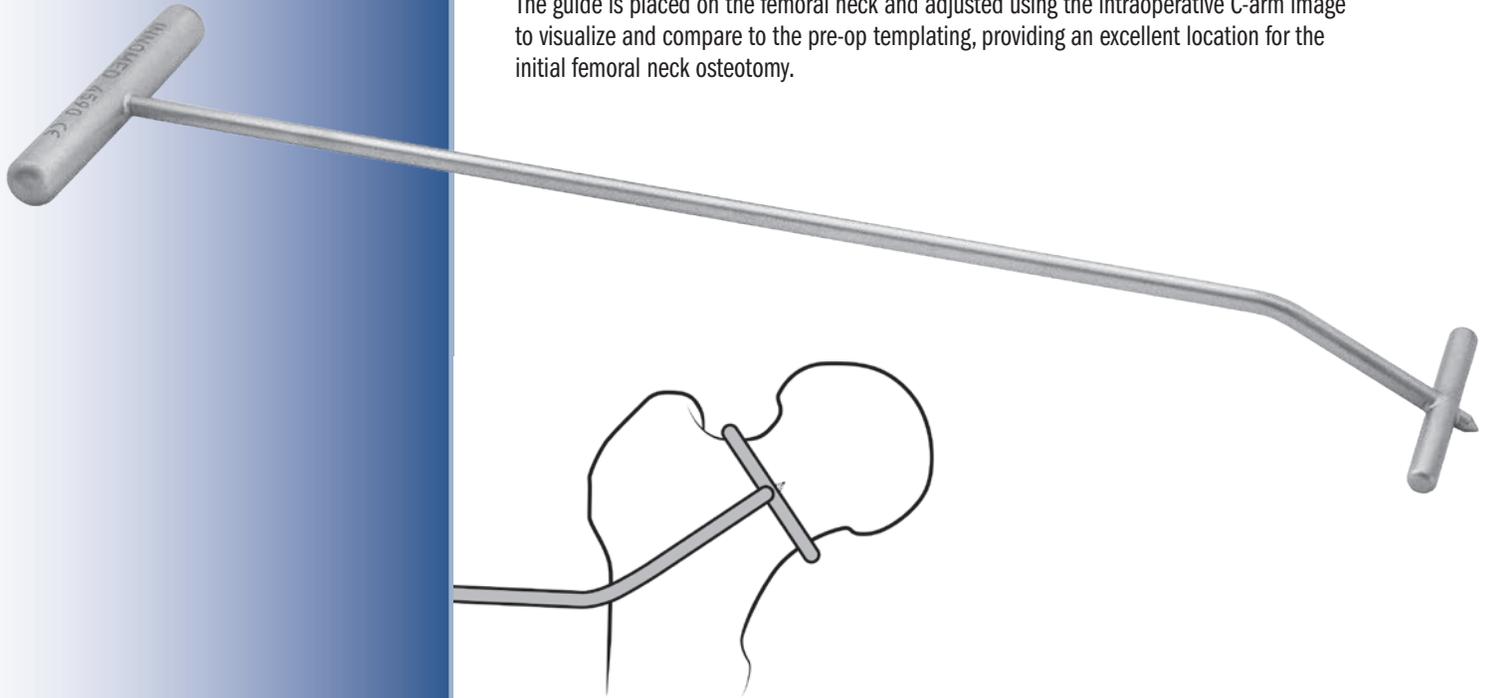

#### PRODUCT NO:

4590

Overall Length: 8.25" (21 cm)  
Handle Length: 1.9" (4.8 cm)  
Cutting Guide Bar Length: 1.22" (3.1 cm)  
End of Bar to Tip Length: 3.5 mm  
Shaft Angle at End: 30°  
Shaft Diameter .125" (3.2 mm)



[www.innomed.net](http://www.innomed.net)

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

# INNOMED



103 Estus Drive, Savannah, GA 31404  
info@innomed.net

912.236.0000 Phone  
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74  
info@innomed-europe.com

1.800.548.2362