

Ludloff/Mau Osteotomy Fixation Clamp

Designed by A. Austin


Used after lateral hallux valgus correction of the metatarsal, the clamp allows for osteotomy fixation and cannulated screw guide wire direction

Clamp fixates the osteotomy to hold the correction, and the 15° slanted cannulated k-wire guide allows the surgeon to place the guide wire for the cannulated screw perpendicular to the osteotomy for final fixation of the osteotomy.

PRODUCT NO:

1812

Overall Length: 5"

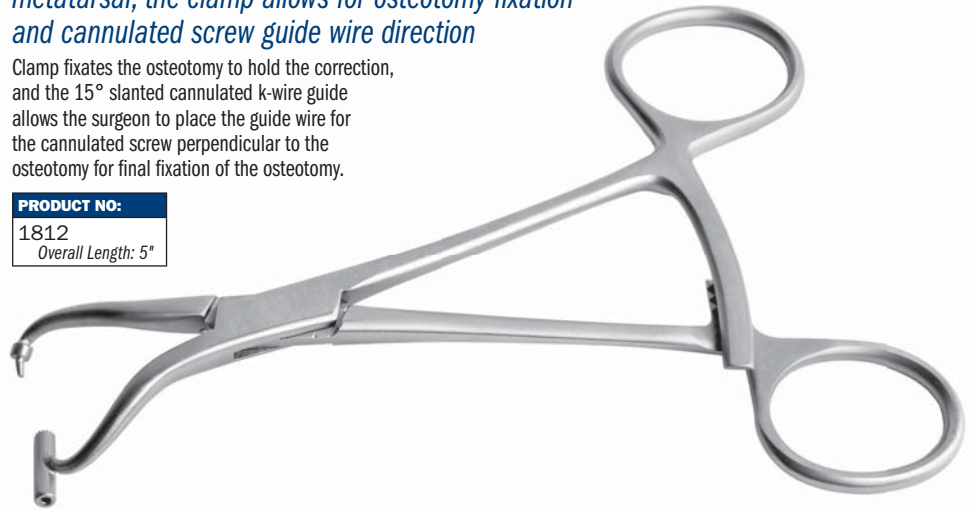

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

© 2011 Innomed, Inc.

Innomed-Europe Tel. +41 41 740 67 73
Fax +41 41 740 67 74

INNOMED

1.800.548.2362