

Meckel Posterior Stabilizing Knee Retractor

Designed by Christopher M. Meckel, MD

Designed to provide enhanced anterior translation of the tibia when doing posterior stabilized total knee replacement

The 15mm deep blade section of the retractor is used to lever the tibia forward (by resting the tip on the posterior tibia and the middle blade section levering off the distal femur) after the box cut has been made in the distal femur.


PRODUCT NO:

4538

Blade Width: 20mm

Blade Depth: 15mm

Overall Length: 10"

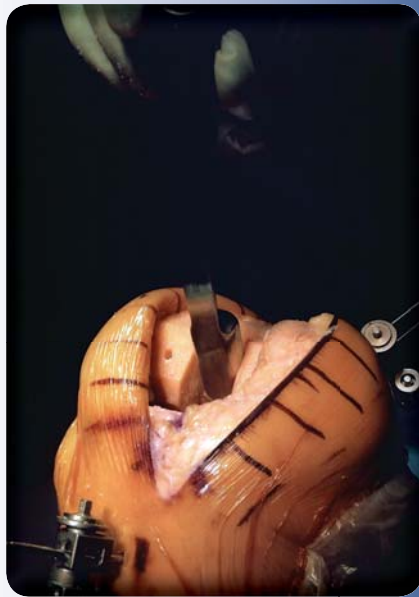

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

© 2011 Innomed, Inc.

Innomed-Europe Tel. +41 41 740 67 73
Fax +41 41 740 67 74



1.800.548.2362