

Screw Extractor Set

Designed to help remove screws with stripped or damaged heads

- ▶ Extractors must be used with drill in reverse.
- ▶ Screw head is reamed with burnishing end, and is then removed with the left turn thread end.
- ▶ Care must be taken to burnish no more than 1/16" (1.5mm) deep, as burnishing too deep can weaken the screw head.

PRODUCT NO:

7250-00 [Set with Case]

7250-01 [Small]
Overall Length: 6"

7250-02 [Medium]
Overall Length: 6"



7250-03 [Large]
Overall Length: 6"



Burnishing Ends



Left Turn Thread Ends

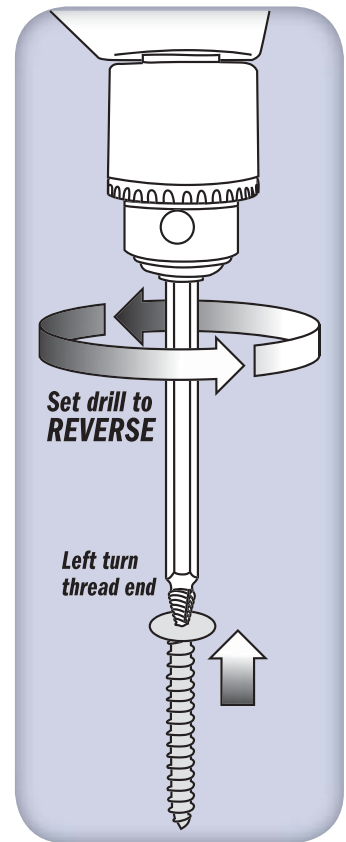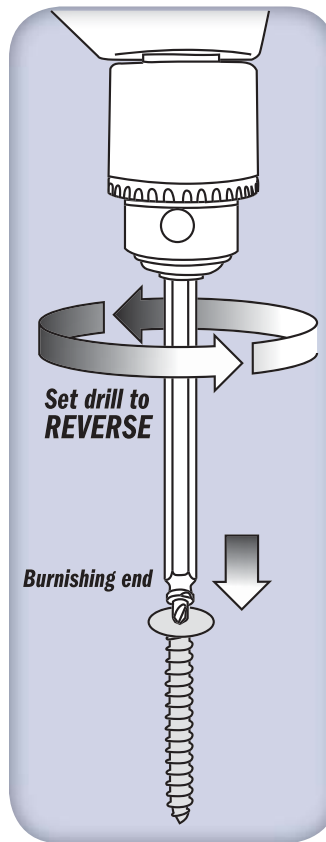

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax



INNOMED

1.800.548.2362