




Gurbani Joint Distractor/Compressor

Designed by Naren G. Gurbani, MD

Versatile joint distractor/compressor provides 360° freedom for arthroscopic or open procedures of foot, ankle, hand, and wrist joints

The surgeon puts the pins in the bone, then slides the holes of the device over the pins and distracts or compresses—the device can be locked in either direction. Especially useful for arthroscopy of subtalar, talo-navicular, calcaneo-cuboid, and wrist joints. The T-wrench helps provide precise, controlled manipulation.

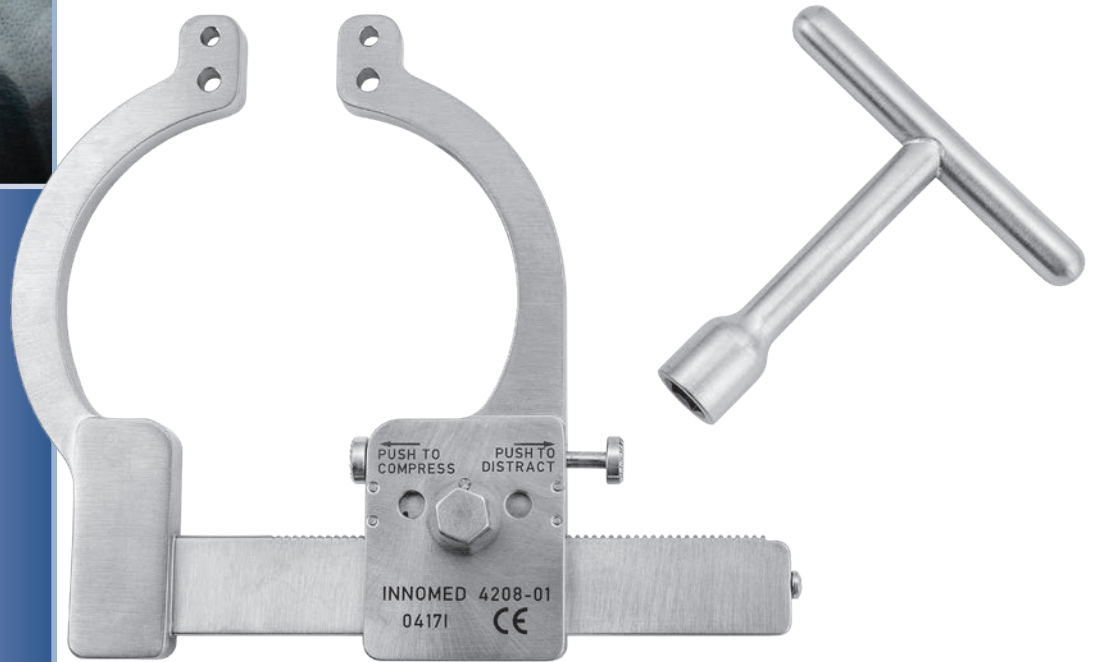

Pin Hole Sizes: .15" (3,5 mm) and .182" (4,5 mm)

| PRODUCT NO'S: |
|--|
| 4208-00 [Set] <i>Includes: Distractor/Compressor, T-Wrench, and Case</i> |
| Available individually: |
| 4208-01 [Distractor/Compressor Only] <i>Dimensions: 6" w x 5" h (15,2 cm x 12,7 cm)</i> <i>Distracts up to: 3" (7,6 cm) / Compresses down to: .5" (1,3 cm)</i> |
| 4208-TW [T-Wrench] <i>Dimensions: 3" w x 3" h (7,6 cm x 7,6 cm)</i> |
| 1025 [Sterilization Case] |



ISO 9001:2008 • ISO 13485:2003

FREE TRIAL ON MOST INSTRUMENTS

Scan to
Launch Our
Website



INNOMED 



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71

1.800.548.2362