



Lateral Condyle Fracture Set

Designed by Carl R. Weinert, MD

Designed for adult and pediatric lateral condyle fractures


Elbow Retractor

Designed for use within the elbow joint to retract the anterior capsule, and provide full exposure of the anterior articular surface for reduction and fixation of displaced lateral condyle fractures. The small blunt tip hooks over the intact medial condyle.



Symmetric Reduction Clamp


Useful to compress T-condylar fractures, and in many other fracture reduction applications. The stops on each end help prevent excessive penetration of metaphyseal and soft bone.

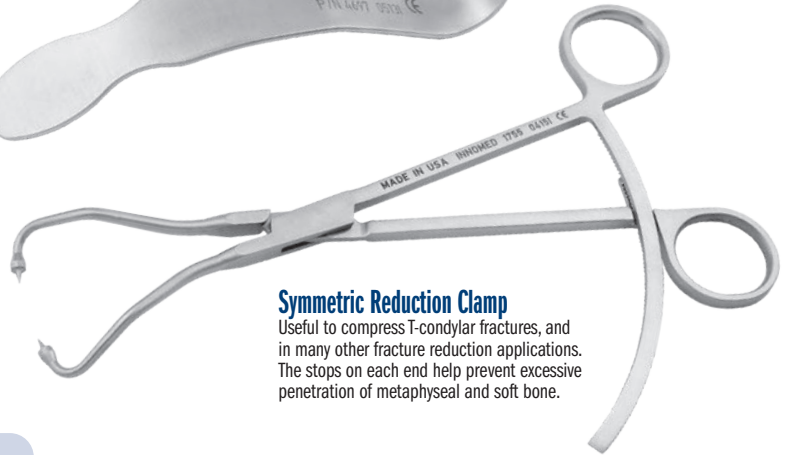

The stops on each clamp end help prevent excessive penetration of metaphyseal and soft bone

Asymmetric Reduction Clamps - Left and Right

Shaped to secure the lateral condyle fragment. The straight tip is placed in the coronoid fossa and the curved tip is used to grasp and compress the lateral condyle fragment. The stops on each end help prevent excessive penetration of metaphyseal and soft bone.

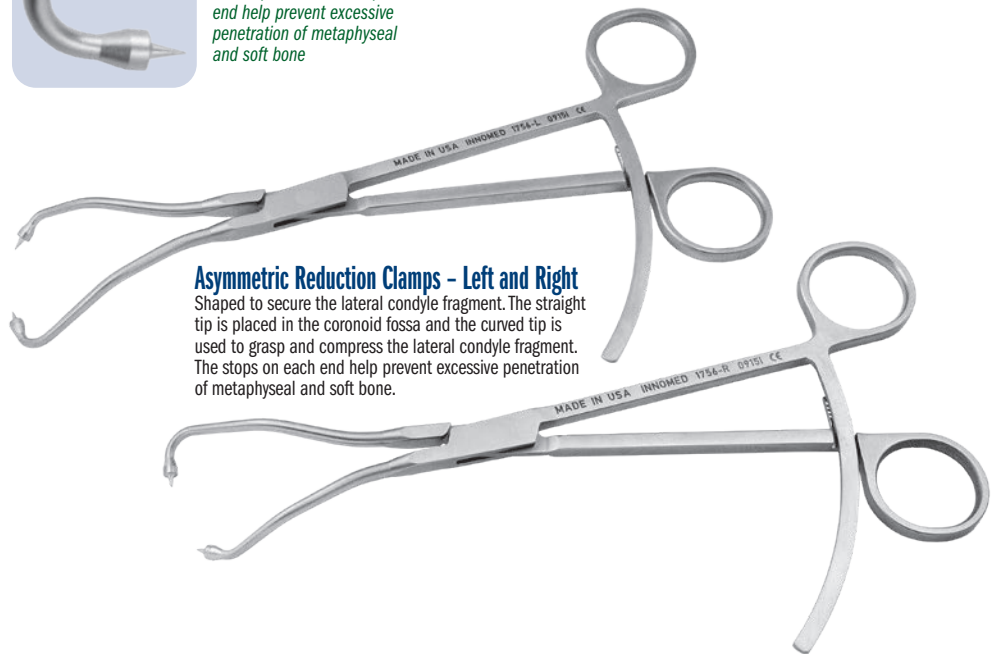

PRODUCT NO'S:

4697-00 [Set with Case]
Set Includes / Available Individually:
1755 [Clamp - Symmetric] Overall Length: 8.5" (21,6 cm) Jaw opens to: 3" (7,6 cm)
1756-L [Clamp - Asymmetric Left] Overall Length: 8.75" (22,2 cm)
1756-R [Clamp - Asymmetric Right] Overall Length: 8.75" (22,2 cm)
4697 [Elbow Retractor] Overall Length: 6.5" (16,5 cm) Blade Width: 1" (2,54 cm)
1015 [Sterilization Case] Dimensions: 11.25" x 7.125" x 3.125" (28,6 cm x 18,1 cm x 7,9 cm)



www.innomed.net

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

INNOMED 



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com

1.800.548.2362

© 2023 Innomed, Inc.